

# AD302L

Lufttrockner

Artikelnummer: 40-1220




## Beschreibung

Der AD302L ist ein praktisches Zuberhörteil für Laboratorien mit externer Luftquelle. Es entfernt Feuchtigkeit, die Ihre Geräte verunreinigen und Resultate beeinflussen kann. Der kompakte, leichte AD302L lässt sich an der Wand anbringen und nimmt so nur wenig Platz weg. Aufgrund der geringen Anzahl beweglicher Teile ist zudem der Wartungsaufwand für das Gerät minimal.



## Generator Details

- **Gasarten:** Null-Luft
- **Integrierter Kompressor:** Nein

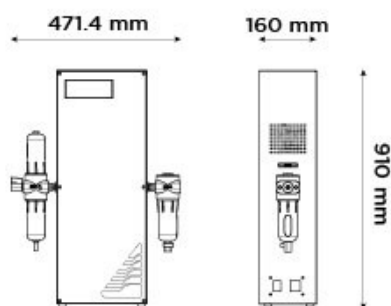

## Hauptmerkmale

- Produziert saubere, trockene Luft und sorgt für eine optimale Versorgung der Gasgeneratoren und -instrumente
- Der Trockner arbeitet bei Bedarf rund um die Uhr
- Der Trockner kann für kleine Labore an der Wand montiert werden.
- Wenige bewegliche Teile bedeuten wenig Wartung und gewährleisten eine lange Lebensdauer des Generators
- Im Betrieb völlig geräuschlos
- Minimale Einrichtung erforderlich
- 12 Monate umfassende Vor-Ort-Garantie

## Technische Spezifikationen

### AD302L - Lufttrockner

Gasart	Null-Luft
Maximaler Gasfluss:	302L/min
Gasanschluss:	1x 1/4" BSP Female
Minimaler Eingangsdruck:	100psi/6,9bar
Maximaler Eingangsdruck:	120psi/8,3bar
Startzeit:	90 mins
Energieverbrauch:	8 - 19.2 W
Stromspannung:	100 - 240 V
Frequenz:	50 / 60 Hz
Strom:	0.8 A
Akkreditierungen:	CE
Größe (HxBxT) mm	910 x 160 x 471 mm
Größe (HxBxT) Zoll	35,8 x 6,3 x 18,5 inches
Gewicht des Generators	20,5kg / 45,1lbs

## [PEAK Protected]<sup>TM</sup>

Peak Scientific gas generators define the benchmark in reliability, convenience and performance in laboratories around the world, and come backed with a 12 month warranty. Beyond this period however you can ensure that your investment continues to be **[Protected]** by our comprehensive generator care cover.

Our world-class aftercare support packages deliver a program of scheduled preventative maintenance whilst giving you the reassurance of instant access to worldwide technical support and priority on-site response in the untimely event of a breakdown.

- For ordering parts visit: [www.peakscientific.com/ordering/](http://www.peakscientific.com/ordering/)
- For service plans visit: [www.peakscientific.com/service/service-plans/](http://www.peakscientific.com/service/service-plans/)

**Peak Scientific UK**  
Tel: +44 (0)141 812 8100  
Fax: +44 (0)141 812 8200

**Peak Scientific North America**  
Tel: +1-800-767-6532  
Fax: +1-978-608-9503

**Peak Scientific China**  
Tel: +86 21 5079 1190  
Fax: +86 21 5079 1191  
For a full list of our worldwide office locations, please visit:  
Web: [www.peakscientific.com](http://www.peakscientific.com)  
Email: [discover@peakscientific.com](mailto:discover@peakscientific.com)

Peak Scientific's Quality Management System conforms to: ISO:9001

