

ZA180A

Stellt Null-Luft als Detektorgas

Artikelnummer: 30-1660 (230v) - UK |
30-1670 (230v) - EU | 30-1680 (110v) -
US




Beschreibung

Der ZA180A stellt Null-Luft als Detektorgas für GC-Anwendungen her und entfernt sämtliche Kohlenwasserstoffe aus einer Druckluftquelle, bevor die Luft in Ihre analytischen Instrumente eintritt. So werden die saubersten Detektorresultate gewährleistet.



Generator Details

- **Gasarten:** Null-Luft
- **Integrierter Kompressor:** Ja

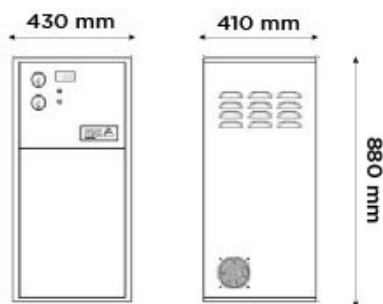

Hauptmerkmale

- Erzeugt keine Luft mit einem Inhalt von $\leq 0,1 \text{ ppm}$
- Führt bei Bedarf rund um die Uhr durch
- Ein Generator kann mehrere Instrumente versorgen
- Wenige bewegliche Teile bedeuten wenig Wartung und gewährleisten eine lange Lebensdauer des Generators
- Minimale Einrichtung erforderlich
- 12 Monate umfassende Vor-Ort-Garantie

Technische Spezifikationen

Der ZA180A stellt Null-Luft als Detektorgas

Gasart	Null-Luft
Maximaler Gasfluss:	18L/min
Kohlenwasserstoffgehalt:	0.1ppm
Maximale Reinheit:	0.1ppm THC
Startzeit:	40 Minutes
Energieverbrauch:	770 Watts
Stromspannung:	110 / 230 ± 10% VAC
Frequenz:	60 Hz
Strom:	2/4 Amps
Maximale Betriebstemperatur:	30°C / 86°F
Akkreditierungen:	CE
Größe (HxBxT) mm	880 x 430 x 410 mm
Größe (HxBxT) Zoll	34,6 x 16,9 x 16,1 inches
Gewicht des Generators	64kg / 140,8lbs

[PEAK Protected]TM

Peak Scientific gas generators define the benchmark in reliability, convenience and performance in laboratories around the world, and come backed with a 12 month warranty. Beyond this period however you can ensure that your investment continues to be **[Protected]** by our comprehensive generator care cover.

Our world-class aftercare support packages deliver a program of scheduled preventative maintenance whilst giving you the reassurance of instant access to worldwide technical support and priority on-site response in the untimely event of a breakdown.

- For ordering parts visit: www.peakscientific.com/ordering/
- For service plans visit: www.peakscientific.com/service/service-plans/

Peak Scientific UK
Tel: +44 (0)141 812 8100
Fax: +44 (0)141 812 8200

Peak Scientific North America
Tel: +1-800-767-6532
Fax: +1-978-608-9503

Peak Scientific China
Tel: +86 21 5079 1190
Fax: +86 21 5079 1191
For a full list of our worldwide office locations, please visit:
Web: www.peakscientific.com
Email: discover@peakscientific.com

Peak Scientific's Quality Management System conforms to: ISO:9001

